FlexiFit™
SERIES

HC 431
FULL FACE MASK

INSTRUCTIONS FOR USE

185 043 239 Rev E

Fisher & Paykel
HEALTHCARE

FISHER & PAYKEL HEALTHCARE

INTERNATIONAL PO Box 14 348, Panmure, Auckland 1134, New Zealand
Tel: +64 9 574 0100
Fax: +64 9 574 0158
Email: info@fphcare.com
Website: www.fphcare.com

USA Tel: 1800 446 3908
or +1 949 470 3900
Fax: +1 949 470 3933

AUSTRALIA Tel: +61 3 9879 5222
Fax: +61 3 9879 5232

AUSTRIA Tel: 0800 29 31 23
Fax: 0800 29 31 22

BENELUX Tel: +31 78 644 0924
Fax: +31 78 644 0914

CHINA Tel: +86 (20) 8363 4286
Fax: +86 (20) 8363 3595

FRANCE Tel: +33 (0)1 64 46 52 01
Fax: +33 (0)1 64 46 52 21

GERMANY Tel: +49 7182 93777 0
Fax: +49 7182 93777 99

INDIA Tel: +91 80 5123 6041
Fax: +91 80 5123 6044

IRISH REPUBLIC Tel: +1800 409 011

ITALY Tel: +39 06 7839 2939
Fax: +39 06 7814 7709

JAPAN Tel: +81 (3) 3661 7205
Fax: +81 (3) 3661 7206

SPAIN Tel: +34 902 013 346
Fax: +34 902 013 379

SWITZERLAND Tel: 0800 83 47 63
Fax: 0800 83 47 54


UK (UK Authorized Representatives) Tel: +44 1628 626 136
Fax: +44 1628 626 146

NTHN IRELAND Tel: +0800 132 189

Fisher & Paykel
HEALTHCARE

FlexiFit™ Series - HC431 Full Face Mask Instructions

The FlexiFit™ HC431 Full Face Mask is intended for multiple patient or single patient adult use by individuals who have been diagnosed by a physician as requiring CPAP or Bi-level ventilator therapy.

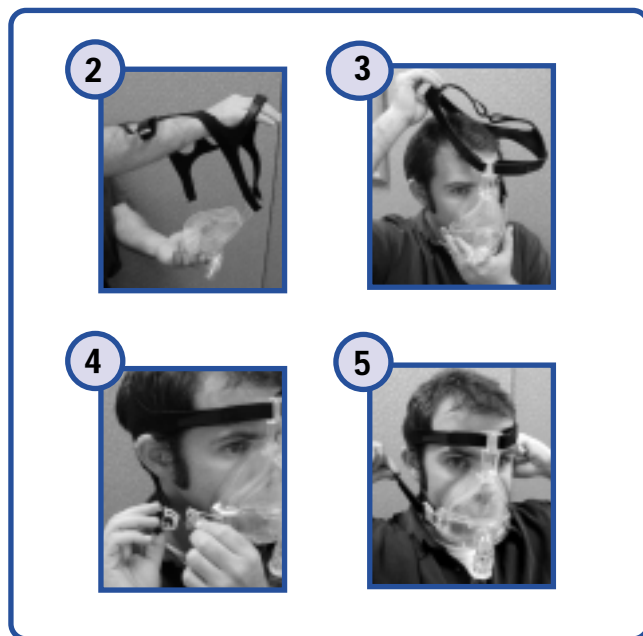FITTING YOUR MASK

Refer to the mask assembly diagram

- Remove plastic cover from Silicone Seal.
- Uncle the Headgear (K) from the Glider™ strap (E) by pulling the Release Cord (I) down. Holding your Mask Base (A) in one hand, spread the Headgear with your free hand just enough to fit comfortably over your head.
- Lightly hold the mask against your face, with the bottom of the seal under your chin, and spread the Headgear over your head with your free hand.
- Re-attach the Headgear Clip (H) to the Glider™ strap.
- Gently adjust the Top Horizontal straps, then the Lower Horizontal straps, then the Top Crown straps. Ensure the top of the seal remains stable on the bridge of the nose, and the bottom of the seal remains fixed under the chin. Ensure that the Glider™ strap is centered.
- Connect the CPAP machine hose running at the prescribed pressure to the mask.
 - If necessary pull the mask slightly away from your face while the straps are still attached. This allows the silicone seal to inflate. Put the mask back on your face.
 - If leaks occur around the upper half of the mask gently tighten the top horizontal straps. It is not necessary for the Forehead Cushion to touch your forehead if it does not.
 - If leaks occur around the lower half of the mask gently tighten the lower horizontal straps. (Avoid over tightening the straps.)
 - If leaks persist, change to the alternate size Silicone Seal, which is found in your FlexiFit™ HC431 Mask box.
- When you need to remove your mask unclip the Headgear from the Glider™ strap by pulling the Release Cord forward or down.

CLEANING YOUR MASK

To ensure continuing safety and comfort for your FlexiFit™ HC431 Mask:

- Before each use:**
- Inspect the mask for deterioration. Do not use if aged or damaged.
 - Inspect the Non-Rebreathing Valve (J). Ensure the silicone flaps occlude the pipe when the CPAP is off.

- After each use:**
- Disassemble the Silicone Seal (C), Elbow (B), Non-Rebreathing Valve (J), and the Headgear (K) from the Mask Base (A).
 - Wash the Silicone Seal, Elbow, Non-Rebreathing Valve and Mask Base in warm soapy water then rinse in fresh water.
 - Leave all parts to dry, away from direct sun light before reassembling.

- Every 7 days:**
- Hand wash the Headgear in lukewarm water using a mild detergent. Allow to air dry out of direct sun light. It is not necessary to remove the Clip (H) from the headgear when washing.

Caution: Do not clean the mask in the dishwasher. Do not clean the mask with alcohol, chlorine or bleach. Do not store the mask in direct sunlight. The above actions may deteriorate or damage the mask and shorten its life. If your FlexiFit™ HC431 weakens or cracks, seek replacement immediately.

ASSEMBLING YOUR MASK

Your FlexiFit™ HC431 Mask is fully assembled and ready to use. However, when you need to clean your mask follow these instructions for reassembly.

- Place Silicone Seal (C) into Mask Base (A).
- Push the Non-Rebreathing Valve (J) onto the Elbow (B), and attach to Mask Base.
- Lay the Headgear (K) flat with the label side down. Position the assembled mask over the Headgear and attach the four straps in to the corresponding slots in the Mask Base and the Glider™ strap (E). This can be done without undoing the Velcro™ Tabs by sliding the Headgear into the slots.

The FlexiFit™ HC431 and accessories are Latex free. Many people find CPAP therapy causes dryness of the nose and throat. Heated humidification can substantially reduce these symptoms. Contact your supplier and ask about Ambient Tracking Plus™.

OXYGEN / PRESSURE PORT CONNECTOR

If pressure readings and/or additional oxygen are required, an Oxygen/Pressure Port Connector is available. (Code 900HC452)

- Connect the large end of the Connector (T) to the Non-Rebreathing Valve (J) on the mask and the small end of the Connector to the Breathing Tube (R).
- Open the Port Caps (P) and place the Pressure and/or Oxygen Tube (S) firmly over the Ports (Q).

(The two ports are identical; either Port may be used for oxygen or pressure measurement. Only remove the Port Cap from ports that are being used.)

- Note: At a fixed flow rate of supplemental oxygen flow, the inhaled oxygen concentration will vary, depending on the pressure settings, patient breathing pattern, mask selection, and the leak rate.
- Warning: If oxygen is used with this CPAP machine, the oxygen flow must be turned off when the CPAP machine is not operating.
- Explanation of Warning: When the CPAP device is not in operation, and the oxygen flow is left on, oxygen delivered into the ventilator tubing may accumulate within the CPAP machine enclosure. Oxygen accumulated in the CPAP machine enclosure will create a risk of fire.
- Warning: Do not smoke while using or in the vicinity of someone using oxygen with this mask.

WARNINGS

- This CPAP mask should be used only with CPAP systems recommended by your physician or respiratory therapist. A mask should not be used unless the CPAP system is turned on and operating properly. The vent slots associated with the mask should never be blocked.
- Explanation of Warning: CPAP systems are intended to be used with special masks with connectors which have vents holes to allow continuous flow of air out of the mask. When the CPAP machine is turned on and functioning properly, new air from the CPAP machine flushes the exhaled air out through the attached mask exhalation port. However, when the CPAP machine is not operating, enough fresh air will not be provided through the mask, and exhaled air may be rebreathed. Rebreathing of exhaled air for longer than several minutes can in some circumstances lead to suffocations.
- At low CPAP pressures the flow through the exhalation port may be inadequate to clear all exhaled gas from the tubing. Some rebreathing may occur.
- Operating pressure range is 3-25 cm H₂O.
- Do not use the mask without the Non-Rebreathing Valve in place. This valve ensures CO₂ rebreathing is kept to a minimum in the event of power failure.
- Do not use the mask if vomiting or feeling nauseous.
- Discontinue using the mask if there is an allergic reaction to any part of the mask. Consult your physician if this occurs.
- If re-using the FlexiFit™ HC431 between different patients, refer to the Multi Patient Cleaning Instructions for the FlexiFit™ HC431 mask.
- Replace the Non-Rebreathing Valve every 12 months.
- Do not use the FlexiFit™ HC431 elbow on other Fisher & Paykel Masks.
- Use the FlexiFit™ HC431 mask only with the headgear provided.

CAUTION

- Use the mask for its intended use only as directed in this booklet.
- US Federal Law restricts this device to sale by, or on order of a physician.
- Operating temperature range is 40-140° F (5-50° C).

CONTRAINDICATIONS

This mask should also not be used by anyone incapable of removing the mask unaided, in the case of vomiting and aspiration, without continuous supervision.

US Patent Number 6 615 834
Other Patents Pending


BIAS FLOW INFORMATION

The FlexiFit™ HC431 has a unique bias diffuser system which vents the expelled CO₂ through two line vents above and below the elbow. It is important that these vents are not blocked by any object. This controlled leak ensures all exhaled CO₂ is expelled from the mask.

Pressure (cm H ₂ O)	3	5	7	9	11	13	15	17	19	21	23	25
Bias Flow (LPM)	23	27	32	37	41	44	48	50	54	57	60	63

RESISTANCE TO FLOW

Pressure drop through mask @ 50 l/min	0.25 cm H ₂ O
Pressure drop through mask @ 100 l/min	0.48 cm H ₂ O

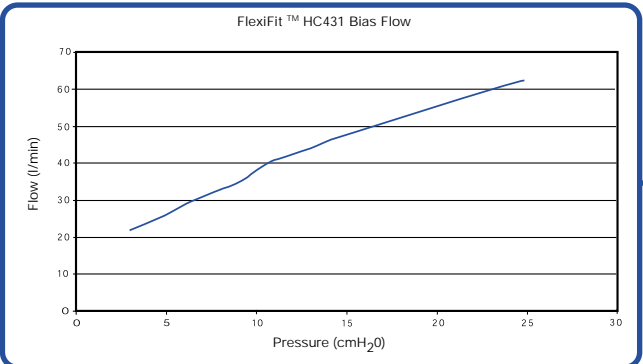FlexiFit™ Series - HC431 Full Face Mask Multi-Patient Cleaning Instructions

Your FlexiFit™ Series HC431 Full Face Mask is provided patient ready. To reuse between patients please follow the guidelines below.


CLEANING / DESOLING

Before disinfecting, disassemble all components from the mask base (The Glider™ strap and Glider Cover do not detach from the Base). Clean all surfaces of mask and components thoroughly with a mild detergent solution and a soft bristled brush. Ensure all visible soil is removed. Rinse thoroughly under running potable water to remove detergent solution.

FAILURE TO PROPERLY CLEAN THIS DEVICE MAY RESULT IN INADEQUATE DISINFECTION

HIGH LEVEL DISINFECTION PARAMETERS

Component	Thermal Disinfection	CIDEX® OPA
Mask Frame	Soak parts in hot water bath for: 30 minutes at 70°C (158°F) or 30 minutes at 75°C (167°F)	Follow manufacturer instructions for high level disinfection. Immerse for 30 minutes at 20°C (68°F). Rinse thoroughly in potable water.
Silicone Seal		
Elbow		
Non-Rebreathing Valve		

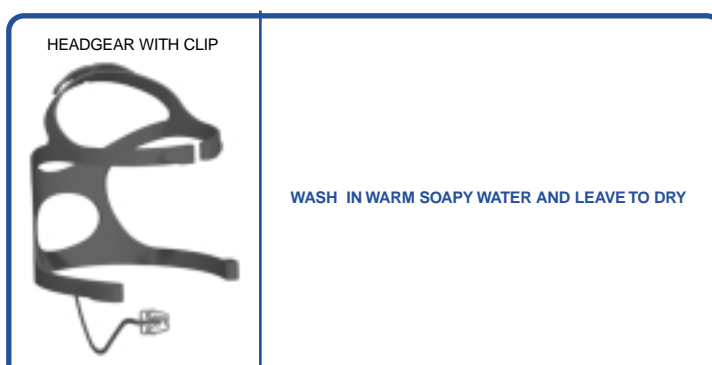WASH IN WARM SOAPY WATER AND LEAVE TO DRY

The FlexiFit™ Series HC431 Mask has been validated to withstand 20 cycles of the above parameters. Do not use if badly damaged. Disclaimer: It is the user's responsibility to qualify any deviations from the recommended methods of reprocessing. Contact your local homecare dealer for additional information.

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